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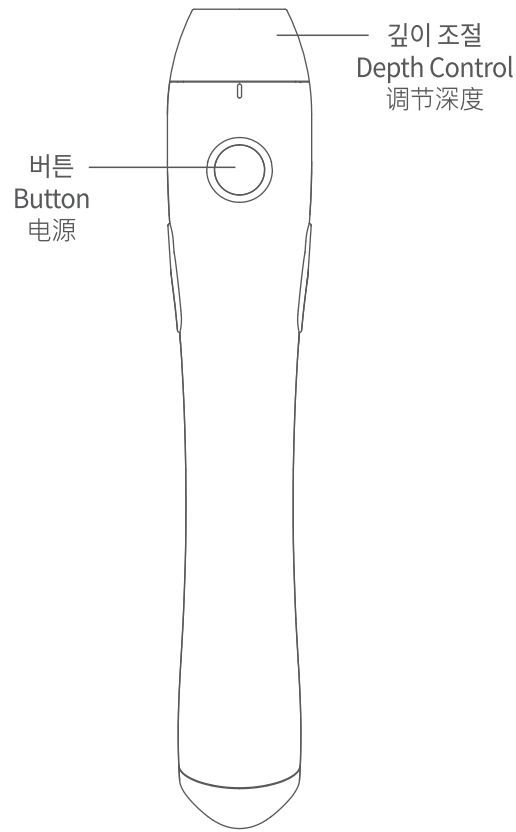
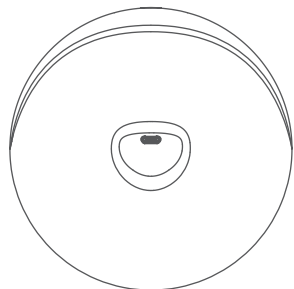
## ■ SPEC

Product name	Skin Beauty Device
Manufactured in	Republic of Korea.
Manufactured by	MCURE Co.,Ltd
Address	4-203, 42-10, Taejunggongdan-gil, Wonju-si, Gangwon-do, Republic of Korea.
Usage duration	2hour, Lv 1 (Continuous use) / Recharge duration : 1h 30m
Product size	Handpiece - W 200mm * D 36mm , Cradle - D 110mm * H 50mm
Purpose of use	This device helps cosmetics to be absorbed with massaging skin.
Battery	lithium polymer (700mAh)
Weight and package unit	1) Weight after packaging : 400g 2) Package unit: 1Set(box)
Efficacy and directions for use	Please refer to the user's manual
Cautions for use	Please refer to the user's manual

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**Product Name** Skin Beauty Device  
**Model Name** AMTS2

Please refer to the S/N on the device

Please refer to the S/N on the device

**Rating**

Input - AC 100~240V 50/60Hz, Output - DC 5V

Please refer to the instruction manual

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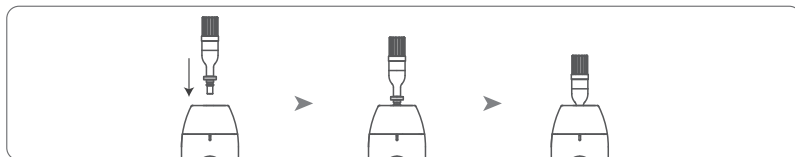
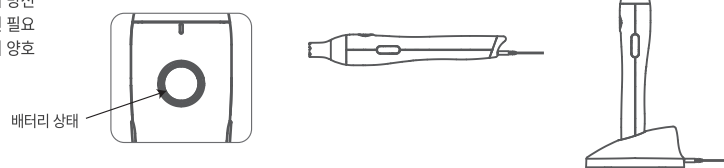
### INSTRUCTIONS BEFORE USE

- 1) Look around if there is any damage on the outside of the package.
- 2) Check if all accessories are contained.
- 3) Check if there is any damage on the cartridge package.
- 4) Use the designated cartridge only.
- 5) Wash your hands and clean the skin.

### INSTRUCTIONS FOR USE

- 1) Please read the user's manual carefully before use.
- 2) Do not disassemble, and/or do not try to remodel this product.
- 3) Make sure cartridge keeps clean when you put into the handpiece.
- 4) Do not touch an adaptor and a power socket with wet hands during charging.
- 5) Do not use this product for children and animals.
- 6) Do not use this product with too much pressure.
- 7) Do not use while device is charging.


### HOW TO KEEP AFTER USE

- 1) Keep it at normal temperature and in the ventilated place.
- 2) Keep out of reach of children.
- 3) This product should be kept with separating the cartridge from the handpiece.

### HOW TO USE


#### ① Put a cartridge into the handpiece.

Put a cartridge into a hand piece, with the cartridge cap covered. Push the cartridge until you can hear the "Tick" sound.
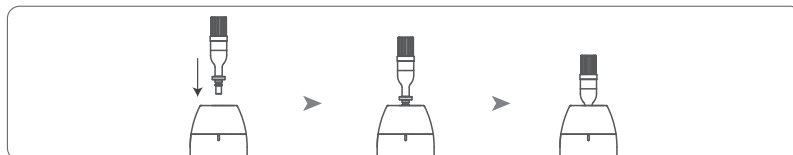
#### ② Put a guide into the handpiece.


Insert the guide into the handpiece, screw it according as picture to the groove around the inside of handpiece, then make sure that guide has locked properly.
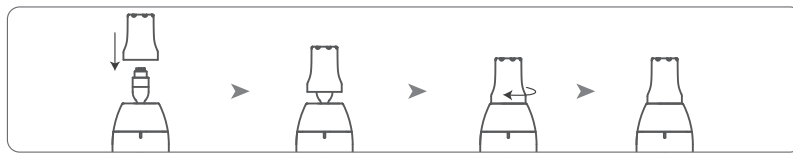
#### ③ Clean the skin.

#### ④ Turn on the device by pressing the button for 2 seconds. Control the operating speed by pressing button.

- Lv 1: 3,000cycles/min
- Lv 2: 4,000cycles/min
- Lv 3: 5,000cycles/min
- Lv 4: 6,000cycles/min

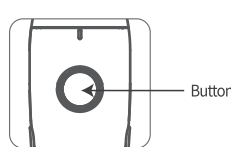
#### ⑤ Once the handpiece equipped with a cartridge and guide massage the skin, cosmetics absorbed through the stimulated area.

#### ⑥ Turn off the device by pressing button for 2 seconds. And then, detach the guide and the cartridge in the written order.

#### ⑦ Charging (Battery state)

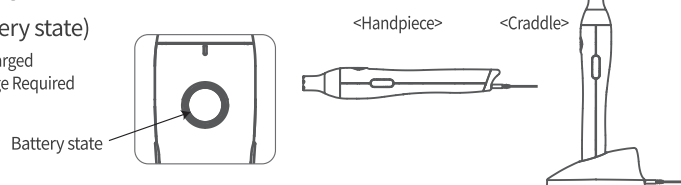- Red: Battery Discharged
- Pink: Battery Charge Required
- Blue: Good

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